


SILICONE SILO BAG - CLEAR CAST

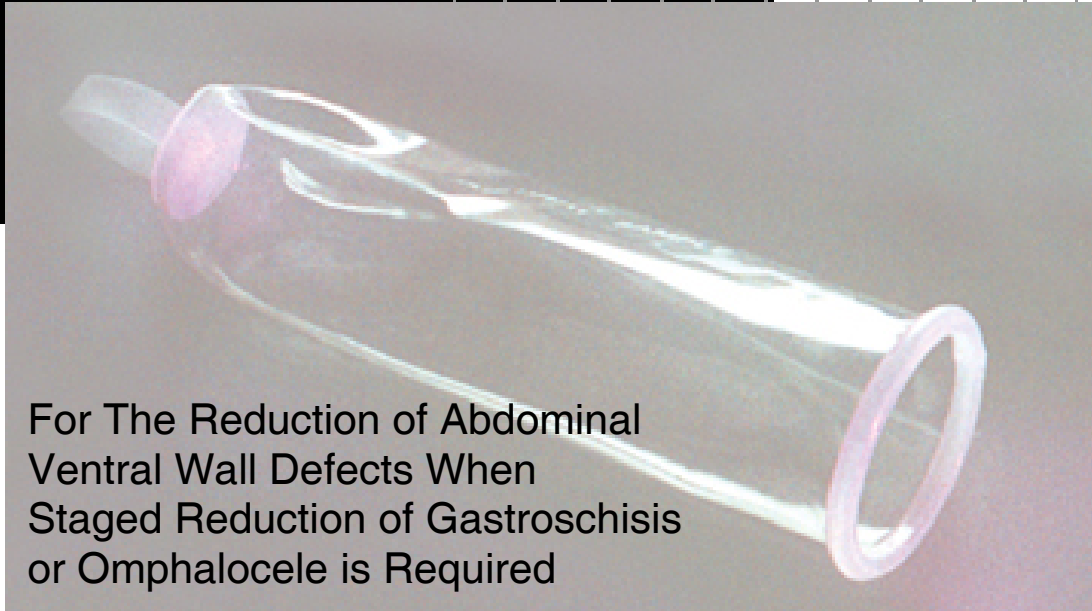

For The Reduction of Abdominal
Ventral Wall Defects When
Staged Reduction of Gastroschisis
or Omphalocele is Required

SPECIFIC DESIGN POINTS:

- **Clear cast silicone elastomer shell** allows for better content visualization.
- **Resilient ring** on open end of bag designed to be placed in the abdomen provides a closed environment.
- **Non-metallic reinforced ring** may be sutured in place when desired.
- **Reinforced silicone loop** located on distal end may be attached to an overhead support for continuous postoperative traction.
- **Six sizes:** 4.0cm, 4.0Wcm, 5.5cm, 6.5cm, 7.5cm, 9.5cm
- **Size imprinted** on each bag.
- **Latex free.**
- **Sizer Rings** to assist in verifying desired silo bag diameter for use.

SILICONE SILO BAG — CLEAR CAST

ORDERING INFORMATION

SUPPLIED STERILE

CATALOG #	DESCRIPTION	DIAMETER	LENGTH
SSB40	Silicone Silo Bag, 4.0	4.0cm	8.0"
SSB40W	Silicone Silo Bag, 4.0 "Wide Body"	4.0cm ring/5.5cm silo	8.0"
SSB55	Silicone Silo Bag, 5.5	5.5cm	8.0"
SSB65	Silicone Silo Bag, 6.5	6.5cm	8.0"
SSB75	Silicone Silo Bag, 7.5	7.5cm	8.0"
SSB95	Silicone Silo Bag, 9.5	9.5cm	8.0"
SSB-SRS	Silicone Silo Bag Sizer Ring Set	All Sizes	n/a


All Silo Bags offered at the same great price of \$200.00/each.

Order for Dr.:	
Patient:	
Surgery Date:	
Catalog #:	Qty:
_____	_____
_____	_____
_____	_____
_____	_____
_____	_____
_____	_____



You asked, and we listened...

Our "wide body" 4.0cm silo design has a 5.5cm silo bag that tapers to the 4.0cm ring, giving you more room to accommodate contents without enlarging the defect.

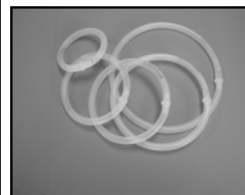

Each Silo Bag is **NOW** packaged with a Sizer Ring to assist you in verifying the silo bag diameter size desired **before** silo bag use!

Limited Warranty, Limitation of Liability and Disclaimer of Other Warranties:

SSP warrants that reasonable care was used in the manufacture and production of this product. Because SSP has no control over the conditions of use, patient selection, surgical procedure, post-surgical stresses, or handling of the device after it leaves our possession, SSP does not warrant either a good effect or against an ill effect following

its use. SSP shall not be responsible for any incidental or consequential loss, damage or expenses directly or indirectly arising from the use of this product. This warranty is in lieu of and excludes all other warranties not expressly set forth herein, whether express or implied by operation of law or otherwise, including, but not limited to any implied warranties of merchantability or fitness for use.



Toll Free 888.878.0811

406.961.0102 · Fax: 406.961.0103 · 1131 North U.S. Hwy 93 · Victor, MT 59875
email: info@ssp-inc.com